



HP-30 PCA

Patient's Comfort Administrator


- ✓ Clear visibility on important information with medical grade large screen.
- ✓ More than 2000 historical events will be recorded for review.
- ✓ Screen lock, prevent unauthorized change in parameter.
- ✓ Drive lock, prevent unauthorized syringe replacement.
- ✓ Labor ward remifentanyl PCA protocol for option.
- ✓ PCEA, PCS function for option.




Simply connect a PCA pack to HP-30 syringe pump, the HP-30 will be turned into a PCA pump.




As a member of HP family, all the advantages of HP family is inherited.

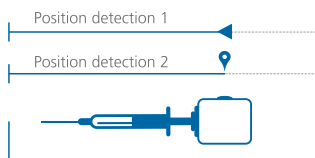

Syringe position detection system.

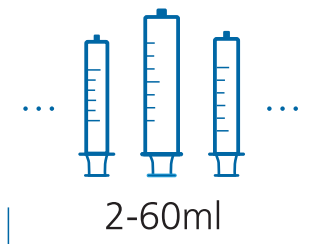

Full range syringe compatibility

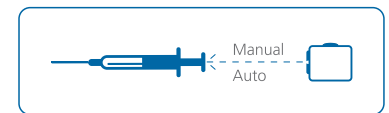

Auto/Manual installation system

Model Type	HP-30 Syringe Pump
Dimensions	258(W)×75(H)×152(D)mm
Weight	1.70kg (including the battery)
Accuracy	≤ 2%
Battery Operating	Time Min. 10h at 5ml/h
Initial Bolus	0 or 0.1-9999.99ml
Basal Rate	0 or 0.1-2000ml/h
Lockout Time	00:00:01-99:59:59 (hh:mm:ss)
PCA Bolus Amount	0.1-9999.99ml
Bolus Rate	0-2000ml/h
PCA Modes	P mode; CP mode; LCP mode

